



Garaventa Stair-Trac



portable wheelchair lift



The Garaventa Stair-Trac is an easy to use, portable wheelchair lift designed to attach under most manual wheelchairs. Stair-Trac's powerful motor allows a small attendant to safely transport a larger passenger up and down stairways. Stair-Trac is suitable for indoor or outdoor use in public buildings and private homes. An affordable and reliable portable access solution.

www.garaventa.ca

STAIR-TRAC - Portable Wheelchair Lift

Specifications

Load Capacity..... 130kg/285lbs
Power Capacity..... 650 stairs
Maximum Speed:
Up Stairs..... 6.5m (21.2ft)/min
Down Stairs 8.8m (28.9ft)/min
Motor 12 VDC
Batteries 2 x 6V - 20AH
Max. Charge Rate.....4 amp
Optimum Stair Angle 35°

Dimensions


Height 922mm/36.3in
Length..... 1468mm/57.8in
Width 636mm/25.0in

Weight

Drive Frame 39.5kg/87lbs
Support Handle 14kg/31lbs
Total Weight 53.5kg/118lbs

Clearances

Wheelchair loading/unloading requires a minimum straight line clearance of approximately 3m/10ft. The diagram below shows typical clearances required on turnback stairways. Dimensions vary with stairway width and operator experience. Consult Garaventa or your local representative for clearance verification.

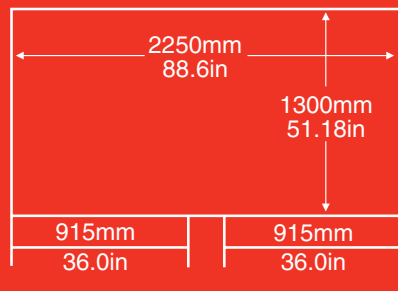

Garaventa Lift

PO Box 1769,
Blaine WA 98231-1769
7505 134A Street, Surrey, BC
Canada V3W 7B3

Phone: 1+(604) 594-0422
Toll Free: 1-800-663-6556
Fax: 1+(604) 594-9915
Email: productinfo@garaventa.ca

Web site: www.garaventa.ca
Global Portal: www.garaventallift.com

©2007 Garaventa. As we are continuously improving our products, specifications outlined in this brochure are subject to change without notice.

15803-E-PB

Features

- Attaches to most standard wheelchairs with push handles (adult/child)
- Powered to carry over 250 lbs. up approximately 30 floors
- Use on virtually any indoor or outdoor stairway
- Stable and self-supporting
- Durable tracks grip stairs securely
- Detachable support handle for compact storage and easy transport
- Quick-release safety belt
- Key switch restricts use to authorized persons
- Rechargeable batteries and charger

Benefits

- Economical compared to stairway lifts and elevators
- It is portable so it can be used on various stairways within a building
- Can be used when electrical power fails or is not available
- Easily stores in a secure area when not in use
- Requires minimal maintenance
- Available immediately
- Can be used for building evacuation

What our customers say

"The Garaventa Stair-Trac has allowed our District to meet its needs at a fraction of the cost of installing an elevator in our building."

- Batavia City Schools, Batavia, NY

"It has proven to be dependable, safe, and functional for our school district."

- Hillcrest Elementary, American Falls, ID



Attach the Stair-Trac to the wheelchair, recline the passenger and move them to the stairway.



Ascend or descend stairs with the press of a button. The passenger sits in a comfortable, upright position.



Stair-Trac will climb virtually any stair tread material.

Authorized Garaventa Representative